

SunTech® *AccuWin Pro v3*[™]

ABPM SOFTWARE

User Guide

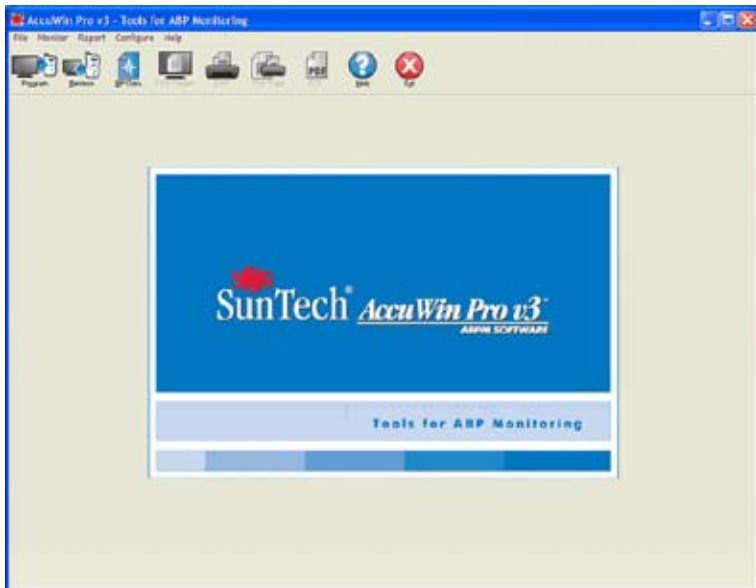


Table of Contents

Introduction to Ambulatory Blood Pressure Monitoring.....	4
<i>AccuWin Pro v3 at a Glance</i>	5
Setting Up the System	6
<i>AccuWin Pro v3 Software Components</i>	6
Hardware Requirements	6
Software Requirements	6
Installing the Software	6
Connecting the monitor to your computer.....	7
Configuring your computer for communication	7
Logging In	8
Entering your ID.....	8
Changing your password.....	8
Changing the current user	8
Using the AccuWin Pro v3 Software	9
<i>AccuWin Pro v3 Basics</i>	9
Toolbar Icons.....	9
Menu Options	10
Programming the Monitor for an ABP Study	12
Retrieving Data from the ABP Monitor	13
Reviewing and Editing an ABP Study	14
Opening a patient file	14
Viewing an ambulatory blood pressure study.....	15
Editing an ambulatory blood pressure study	16
Entering Patient Diary Comments	17
Changing Awake/Asleep Times	18
Setting BP Threshold Limits	18
Setting Pediatric Threshold Limits.....	19
Defining Time-Slice Periods.....	20
Editing Patient Info and Physician Info.....	21
Viewing Statistics.....	23
Viewing Time-Slice Statistics.....	24
Viewing Hourly Averages	24
Viewing the Interpretive Report Summary.....	25
Comparing two studies.....	26
Creating Reports	27
Configuring and customizing the report	27
Report - Formats tab	28
Report - Details tab.....	29
Previewing the report.....	30
Printing the report	30
Saving the report as a PDF	30

Managing Patient Files 31

- Open..... 31
- Copy 31
- Export 31
- Email 32
- Delete..... 32
- Converting patient data files 33

Customizing and Configuring AccuWin Pro v3 34

- Hardware Settings..... 34
- Data Settings..... 35
- Color Settings..... 36
- Diary Comment Settings..... 37
- Email Settings..... 37
- Report Settings..... 38
- Export Settings..... 39

Overview of AccuWin Pro v3 Report Symbols 40

- Report Symbols..... 40

Overview of Event, Quality, and Test Codes 41

Administrative Tools 43

- Maintaining login options..... 43
- Adding user accounts..... 45
- Changing user level 45
- Deleting user accounts 46
- Resetting user passwords 46

AccuWin Pro v3 Technical Assistance..... 47

Safety 48

Index 50

Introduction to Ambulatory Blood Pressure Monitoring

Ambulatory blood pressure monitoring is an accepted clinical tool for collecting multiple blood pressure measurements to better assist clinicians with diagnosis and management of hypertension by providing: blood pressure variability, estimation of true blood pressure, overnight changes in blood pressure, and morning surge in blood pressure.¹ In-clinic or home blood pressure measurements are not representative of blood pressure readings throughout the 24-hour day. Several studies have shown that ambulatory blood pressure monitoring, when compared to clinic or home blood pressure measurements, is superior in predicting target organ damage, morbid events, or cardiovascular risk.^{2, 3, 4}

The data obtained from ambulatory blood pressure monitors is highly accurate and useful for managing a wide variety of hypertensive situations including:

- White-coat hypertension
- Resistant hypertension
- Masked hypertension
- Childhood hypertension
- Efficacy of anti-hypertensive drug therapy on a 24-hour basis
- Nocturnal hypertension
- Episodic hypertension and/or anxiety disorders
- Hypotensive symptoms
- Changes in diet and daily routine designed to reduce hypertension
- Hypertension in pregnancy

1. Pickering, T.G., Shimbo, D., & Haas, D. (2006). Ambulatory Blood-Pressure Monitoring. *New England Journal of Medicine*, 354(22), 2368 – 2374.

2. Pickering, T.G., Shimbo, D., & Haas, D. (2006). Ambulatory Blood-Pressure Monitoring. *New England Journal of Medicine*, 354(22), 2368 – 2374.

3. Marchiando, R.J. & Elston, M.P. (2003). Automated Ambulatory Blood Pressure Monitoring: Clinical Utility in the Family Practice Setting. *American Family Physician*, 67(11), 2343 – 2350.

4. White, W.B. (1999). Ambulatory blood pressure as a predictor of target organ disease and outcomes in the hypertensive patient. *Blood Pressure Monitoring*, 4(3), 181 – 184.

AccuWin Pro v3 at a Glance

AccuWin Pro v3 is a simple, powerful application that allows maximum flexibility in the configuration, analysis, interpretation, and reporting of ABPM studies. *AccuWin Pro v3* offers the first automated interpretive summary of ABP data with analysis of the white coat effect, advanced online Help system, multiple report configurations, and integrated e-mail capability.

The program's easy-to-use, menu-driven design enables you to analyze ABP data in a number of ways. The interpretive report summary feature assists with evaluating ABP levels for normality. The single page report now provides the results and patient information needed on a single sheet. The hourly averages and study comparison features better display BP trends to assist in patient treatment. Statistical calculations, including pulse pressure and pressure rate product, provide insight into patient activity levels and outcomes.

With *AccuWin Pro v3*, you can design, print or create PDF reports using either the default or your own custom formats. When creating your own report formats, you can include any combination of the following information:

- ABP Summary including interpretive analysis, graph, and statistical results
- Patient information details
- Complete statistical details
- Detailed record of collected ABP readings
- Hourly averages of collected ABP readings with additional statistical analysis
- Study comparison
- Histograms of ABP data and statistics
- Review of programmed ABP monitor configuration

A Note about HIPAA

The regulations set forth by the Health Insurance Portability and Accountability Act of 1996 (HIPAA) cover a number of topics, two of which present possible compliancy issues for health care providers using a software application such as *AccuWin Pro v3*: **privacy** and **security**.

Password-controlled user access and encrypted patient information capabilities of *AccuWin Pro v3* assist health care providers in maintaining a HIPAA compliant environment. Since HIPAA compliancy is ultimately the responsibility of the provider, please be aware that report printouts, report PDF files, and exported data contain unencrypted patient information and should be handled appropriately.

Setting Up the System

AccuWin Pro v3 Software Components

The *AccuWin Pro v3* Software components include:

- *AccuWin Pro v3* User's Guide
- *AccuWin Pro v3* CD
- Oscar 2™ PC USB interface cable

Hardware Requirements

- Pentium-based or equivalent PC-compatible computer with CD drive
- SVGA or compatible display adapter and monitor (1024 x 768 pixels minimum recommended resolution)
- One USB serial port
- 32 MB RAM and minimum 20 MB hard drive recommended

Software Requirements

- Microsoft Windows version 98 or later, excluding Windows ME (see page 47 for instructions on Windows 2000)
- Microsoft Internet Explorer 3 with Microsoft HTML Help

NOTE: *AccuWin Pro* does not support a server license. Multiple users simultaneously operating from the same copy of the *AccuWin Pro* software, either located on a single machine or server, is not supported.

Installing the Software

NOTE: If a previous version of *AccuWin Pro v3* is already installed, we recommend backing up any previously collected patient files before upgrading.

Place the installation CD in the CD drive located on your computer. If CD autoplay is enabled on your computer, follow the instructions that appear on your screen.

If autoplay is not enabled, follow these steps:

1. Open Windows Explorer or Windows NT Explorer (Press the Windows "Start" button and find "My Computer" or "Computer").
2. Click on the **CD drive**.
3. Double-click the **AUTORUN.EXE** file.
4. Follow the instructions on the screen.

NOTE: The USB cable should not be connected to the computer before *AccuWin Pro* is installed.

Connecting the monitor to your computer

The first step in communicating with the device is to connect it to your computer. If using the RS-232 serial cable, make sure that the power to the PC is off to avoid electrical shock.

1. Connect the PC interface cable to the connection at the bottom of the ABP monitor (Fig. 1).
2. Connect the other end of the PC interface cable to a USB port on the back of your computer (Fig. 2).
3. The first time you connect the USB cable to your computer, you may see an "Installing Drivers" message. Wait until the drivers are installed before programming the *Oscar*.
4. If using the serial port PC interface cable, connect the cable to the 9-pin serial port on the back of your computer.



Fig. 1



Fig. 2

NOTE: The cable can be left connected when the PC is off.

Configuring your computer for communication

Installing *AccuWin Pro v3* will load the driver(s) for the USB cable. Once the cable is connected to the PC, *AccuWin Pro v3* will recognize the cable and auto-select it as the connection to the monitor (Fig. 3).

If using a serial port instead of a USB port, the port will need to be set correctly within *AccuWin Pro v3*.

1. From the **Monitor** menu, select **Configure Port**. The hardware settings tab of the configuration window is displayed (Fig. 3).
2. Select the **serial port** to which the PC cable is connected.
3. Select the **ABP device** corresponding to the monitor that you have.
4. To verify that the settings are correct, click **Test Port**. If communication is successful, a confirmation message will be displayed quickly.
5. Select **OK**.

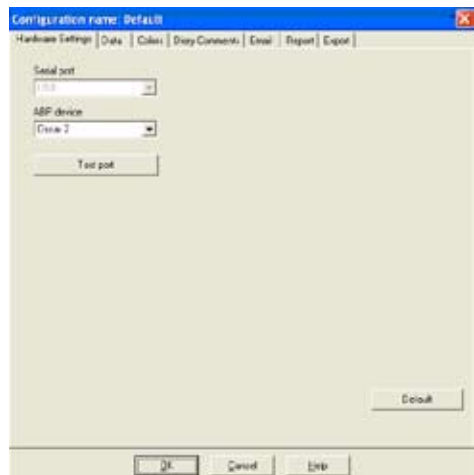


Fig. 3: Hardware Settings

Logging In

AccuWin Pro v3 protects sensitive patient information by letting the system administrator enable a login process that requires all users to log in before accessing any information in *AccuWin Pro v3*. By default, log in is not required unless the system administrator enables this. To enable this feature, Choose **Configure > Admin Tools**.

Entering your ID

If your *AccuWin Pro v3* administrator enables password control, you must enter a valid user name and password when you open *AccuWin Pro v3* or change the current user.

NOTE: The administrator can program *AccuWin Pro v3* to automatically log a user off after a specified time of inactivity elapses. If you are logged off, you must log back in again.

To log in:

1. Enter a valid user name.
2. Enter a valid password.
3. Click OK.

Changing your password

If you know your current password, you can change it. If you are logging in for the first time, you must change the default password assigned to you. You must be logged in under your name before you can change your password.

1. From the Configure menu, choose Change Password.
2. Enter your current password.
3. Enter your new password. It can be alphanumeric and it must be at least 6 characters.
4. Re-enter your new password to confirm it.
5. Click OK.

NOTE: If you forget your password, you must ask the administrator to reset it. However, if you know your current password, you can change it at any time.

Changing the current user

1. From the Configure menu, select Change User.
2. Enter a valid user name.
3. Enter a valid password.
4. Click OK.

Using the *AccuWin Pro v3* Software

AccuWin Pro v3 gives you the flexibility to program the ABP monitor the way you want and the simplicity to collect and retrieve the important test data.

AccuWin Pro v3 Basics

- The toolbar provides easy access to frequently used functions.
- The menu bar allows you to access all functions of the software.

Toolbar Icons

- **Program:** Programs the monitor for an ABP study.
- **Retrieve:** Retrieves an ABP study from the monitor.
- **BP Data:** Opens an ABP study.
- **Print Preview:** Previews a report.
- **Print:** Prints a report.
- **Print Page:** Prints the report page of the displayed data.
- **PDF:** Saves the report as an Adobe Portable Document Format (PDF).
- **Help:** Displays the online help.
- **Exit:** Closes *AccuWin Pro v3*.



Fig. 4: Opening window

Menu Options

AccuWin Pro v3 software includes the following menu options:

File:

- **Open:** Opens a list of patient data files from which you can select.
- **Close:** Closes the currently displayed patient data file.
- **Copy:** Displays a list of patient data files. Select the file you want to copy.
- **Export:** Lets you export *AccuWin Pro v3* data files to files in either ASCII or GDT format. You can choose to export the currently displayed file or choose one from the list of patient files.
- **Email:** Creates a new email message, attaches selected patient data files, and uses Simple Mail Transfer Protocol (SMTP) to send the email to its destination.
- **Delete:** Deletes one or more patient data files. Each patient file contains ALL the BP studies for that patient; therefore, ALL BP studies for that patient will be erased when a patient file is deleted.
- **Exit:** Closes *AccuWin Pro v3*.

NOTE: There are several short-cut keys programmed into *AccuWin Pro v3* that can be identified by a letter underlined. The function can be executed by pressing the "Alt" key and the underlined letter. For example, open *AccuWin Pro v3* help by pressing "Alt" and "h."

Monitor:

- **Program study:** Prepares the monitor for a patient study by letting you set the parameters for your patient.
- **Retrieve data:** Retrieves patient data collected on the ABP monitor.
- **Configure port:** Lets you configure your computer to communicate with your ABP monitor. You can also test the settings that you selected to make sure it is operational.
- **View monitor configuration:** Displays the configuration settings for the monitor used in the displayed study.

Report:

- **Configure:** Lets you choose report formats and customize report settings.
- **Preview:** Displays the report with the currently configured format as it will print.
- **Print:** Prints the report associated with the currently configured format.
- **Print Page:** Prints the page of the report that corresponds with the selected tab in the *AccuWin Pro v3* main window.
- **Create PDF:** Saves the report as an Adobe Acrobat Portable Document Format (PDF). This format captures all font and layout information of the report. It is a format that you can use to send to others for review.

Configure:

- **Preferences:** Lets you configure the way you collect and process BP data. You can change hardware and data settings, customize colors in the graphs and charts, customize a comments list from which physicians can select to record their comments about patients, configure email settings, configure report settings, and set data export options.

If you manage the data for a large medical practice or a service bureau conducting clinical trials for multiple clients, you can generate several configuration profiles—for example, one for each physician in the practice or each service bureau client.

- **Admin Tools:** Lets you create new user accounts, change access levels (administrative or user access), reset user passwords, enable the login process, and enable automatic logging off after a specified time elapses.
- **Change password:** (Appears only if the login process is enabled by the *AccuWin Pro v3* administrator.) Lets you change your own password. You must know your current password to do this.
- **Change user:** (Appears only if the login process is enabled by the *AccuWin Pro v3* administrator.) Lets you change the currently logged in user.
- **Diagnostics:** Lists diagnostic information for *SunTech* technical support.

Help:

- **Help topics:** Displays the online Help.
- **SunTech Web site:** Launches *SunTech's* Web site.
- **About:** Displays the *AccuWin Pro v3* version and your computer's available memory and platform.

Programming the Monitor for an ABP Study

To prepare the monitor for an ABP study, simply fill out an on-screen form to set the parameters for your patient to be programmed into the monitor.

To program the monitor:

1. From the **Monitor** menu, select **Program study**. Alternatively, click the **Program** icon on the toolbar.
2. Enter the settings in the form (Fig. 5).
3. Click **OK**. An indicator bar shows the progress as the data is transferred to the monitor and disappears when programming is successfully completed.

Fig. 5: Program Monitor

The test parameters can be adjusted as follows:

- Patient name and ID:** For reporting and referencing data.
- Start study in 5 minutes:** Check denotes that the study will start automatically after programming; unchecked denotes that the study must be manually initiated with the Start/Stop button.
- Time zone difference:** Adjust the monitor's clock to the time zone that the patient is in relative to your time zone.
- Max Pressure:** 160 to 280 mmHg. Suggested setting is 30 mmHg above the highest expected systolic BP.
- NOTE:** The ABP monitor will not inflate to Max Pressure with each reading; it inflates to 30 mmHg above the previous systolic reading.
- Keypad:** **Enabled** allows the patient to start readings.
- Display:** **Enabled** allows the patient to view the results immediately after a measurement.
- NOTE:** Keypad and Display are always enabled for the first 30 minutes of study.
- Intervals:** Set intervals between programmed readings. Select **Fixed** to set the intervals to exact times. Select **Standard** for +/- 5 minutes around the selected intervals.
- NOTE:** 5 and 10 minute intervals are always fixed.
- Time Periods:** Up to 3 allowed.
- Time Intervals:** 5, 10, 15, 20, 30, 45, 60, 90 and 120-minute intervals between readings.

Retrieving Data from the ABP Monitor

Before retrieving data from the monitor, ensure that the monitor is disconnected from the patient before connecting the ABP monitor to any piece of hardware (such as your computer).

To retrieve the data:

1. Connect the ABP monitor to the computer.
2. From the **Monitor** menu, select **Retrieve Data**. Alternatively, click Retrieve on the toolbar. Communication starts automatically. The dialog box on the screen shows the progress as the data is transferred.
3. After completion, a dialog box asks whether this is a new patient. If you click **No**, you will be presented with a list of patient files already present on your hard drive (Fig 6). Select the patient file to save the retrieved data.
4. If you click **Yes**, you will be presented with the **Patient information** dialog box (Fig 7). Enter data in the fields. If the Enable prompt for file name is set in the data configuration you will be asked to enter the file name.



Fig. 6: Retrieving data, choosing patient

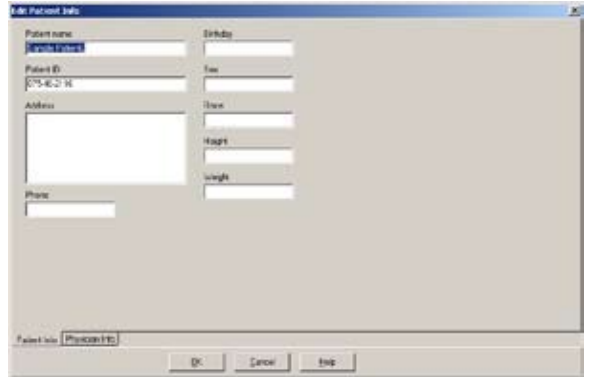


Fig. 7: Retrieving data, patient information

NOTE: If using a "." within the filename, any characters after the final "." will be replaced by ".awp." For example, a filename of john.smith.2008 will be saved as john.smith.awp. To save the file as intended, enter an additional "." at the end of the filename. For example, john.smith.2008. will be saved as john.smith.2008.awp.

5. Click the **OK** button to save the data. The BP file will now open automatically.

The patient's BP data is now saved as a file on your computer. The ABP Data tab displays the retrieved data. The patient's name, ID#, test date, and file name and its location are displayed in the status line located at the bottom of the display area.

Reviewing and Editing an ABP Study

Opening a patient file

1. *AccuWin Pro v3* will automatically open a patient file directly after the data is retrieved.
2. From the File menu, select Open. Alternatively, from the toolbar, click the BP Data button.

NOTE: If you are opening patient files not converted to the current version of *AccuWin Pro v3*, a message appears requesting that you convert previously saved patient files to new encrypted files to enhance security. See *Converting patient data files*.

3. The Open APB Study/Choose Patient dialog box appears, allowing you to select the patient file to open. (An asterisk next to the file name indicates patient files that have not been converted.) Select the patient whose BP data you wish to retrieve by clicking on the patient's name.
4. Select OK.
5. If the patient has more than one ABP study, a second dialog box, Choose Date, appears. Select a date and select OK. The patient's data appears in the display area. You can now edit the BP studies or print a report.

Each study in a patient file contains data that is displayed in the views accessed by the tabs at the bottom of the application. You can edit data on the following tabs: ABP Data, Patient Info, Physician Info and Summary. The status bar identifies the patient name, ID, study date, and file name of the displayed file.

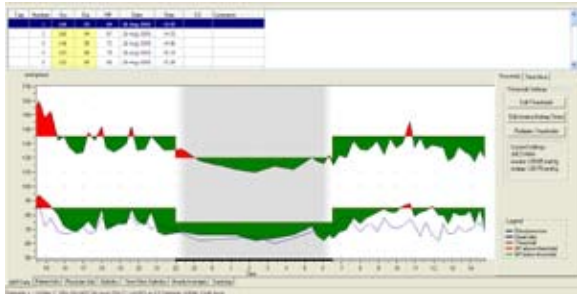


Fig. 8: ABP Data Graph, Threshold tab

Tabs and their contents

- **ABP Data:** ABP study data from the monitor.
- **Patient Info:** Patient name, ID, contact information, and physical description.
- **Physician Info:** Patient history, clinical information, and interpretation.
- **Statistics:** Statistical analysis of the edited ABP study.
- **Time-Slice Statistics:** Statistical analysis of specific user-defined time windows.
- **Hourly Averages:** Analysis of edited ABP study data including study comparison.
- **Summary:** Interpretative summary settings and results for current study.

Viewing an ambulatory blood pressure study

In the upper portion of the ABP Data tab, a table displays each of the BP measurements collected during the test period (Fig. 8). The table shows highlighted systolic and diastolic BP readings, heart rate, and any event, quality, or test codes associated with the measurement.

Below the table is a graphical representation of the BP data. The scale of the vertical axis represents blood pressure in mmHg and heart rate in beats per minute (bpm). The horizontal axis displays the time in clock hours. Clicking on any point in the graph highlights the corresponding measurement row in the table.

There are two tabs to the right of the graph: Threshold and Time-Slice. The Threshold tab displays the legend for the ABP Data graph and allows the user to change the Threshold settings on the graph. The Time-Slice tab allows users to define specific time windows in which to perform a statistical analysis.

Shading on the graph indicates the period of the study. The default colors (which can be changed using the **Configuration > Colors** tab) are the following:

- **Light blue shading:** Indicates the white coat period, which is the first hour of the study (if enabled through the Summary tab)
- **Light gray shading:** Indicates the asleep period of the study.

To zoom in on a portion of the graph:

1. Place the mouse over either axis of the graph (Fig. 8) until it turns into a double arrow.
2. For the vertical (horizontal) adjustment, drag the double arrow down (right) to zoom in or up (left) to zoom out. If zooming in causes part of the graph to be hidden, the right-hand or bottom border turns into a scroll bar that you can use to adjust the viewing window.

To adjust the size of the data table:

1. Position the mouse pointer over the bottom boundary of the table.
2. When the pointer turns to a double arrow, drag it up or down to resize the table and graph displayed.

Editing an ambulatory blood pressure study

In the ABP Data view, you can edit an ABP study for accuracy and context. When the data is retrieved by *AccuWin Pro v3*, all readings in an ABP study that have suspect codes will be “tagged” in the first column of the table, either with an asterisk (*) or an exclamation point (!). Tagged records are omitted from the analysis of the ABP study displayed in the Statistics and Hourly Averages views, but these readings can be printed in the report (see Omitted BP Data Page).

The monitor tags data with an asterisk based on criteria used to determine the validity of the data. The exclamation point (!) is a permanent tag and is used when an accurate reading can not be determined by the ABP monitor. This tag can not be changed and associated data will not be used in data analysis. The asterisk (*) is a tag that you can edit. Data can be tagged or un-tagged based on a number of factors, including patient history, patient diary information, or other factors.

There are two methods to remove or insert an asterisk (*) tag on a measurement.

Using the table:

1. Highlight the reading you want to tag (omit) or un-tag (include) by selecting it.
2. Select the first, or left-most, cell under the column, labeled “Tag”

From the graph:

1. Use the cursor on the graph to select the reading you want to tag (omit) or un-tag (include). Consequently, this action highlights the reading in the table.
2. Right-click the mouse which invokes the context menu.
3. Select Toggle Tag from the menu.

NOTE: All edits to the ABP study data are saved immediately.

Entering Patient Diary Comments

Use the Comments column, the right-most column in the table, to keep track of patient activity during a BP reading. While activity is only one of the many factors that can affect blood pressure, it can be helpful in understanding a BP reading within the context of the study.

To simplify this task for the patient, a patient diary checklist is included with your ABPM system. The same entries are available in the list accessed by the cells under the Comments column.

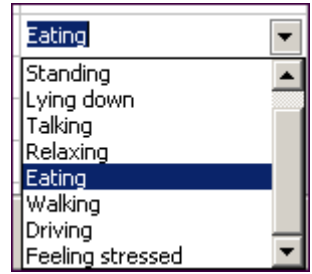


Fig. 9: Comments from table

To enter a patient comment from the table (Fig 9):

1. Highlight the reading to which you want to add a comment, and select the corresponding cell in the Comments column.
2. Press the left mouse button to invoke a pull-down menu.
3. Enter a comment or select a comment from the list.
4. Press Enter.

To enter a comment from the graph:

1. Using the mouse pointer, select the reading in the ABP data graph to which you want to add a comment by pressing the left mouse button. This causes the corresponding row in the table to be highlighted
2. Press the right mouse button, which invokes the context menu.
3. Select Set comment from the menu.
4. Add comment by typing or selecting from list, and select OK to add.

NOTE: All edits to the ABP study data are saved immediately.

You can add or delete entries appearing in the comments list by going to the Diary Comments section of the Configuration window. (For instructions on configuring Diary Comments, see *Customizing and Configuring AccuWin Pro v3.*)

Changing Awake/Asleep Times

Along with BP thresholds, asleep times are used to calculate BP loads and asleep dip percentage. The time is defined by the parameters set when preparing the monitor for a study. However, you can change these times if they have not been pre-defined or are inaccurate.

1. Display an ABP study (Fig. 8).
Choose the Threshold tab.
2. Click Edit Awake/Asleep Times to change the patient's asleep and awake times.
3. In the Times tab (Fig. 10), reset the patient's asleep and awake times to the nearest half-hour by moving the appropriate sliders right or left. The new times are displayed on the clock face, with the black portion representing time asleep and the white portion, time awake. The clock is displayed in military time, i.e., 12:00 is noon, 18:00 is 6 p.m., and 00:00 (or 24:00) is midnight.

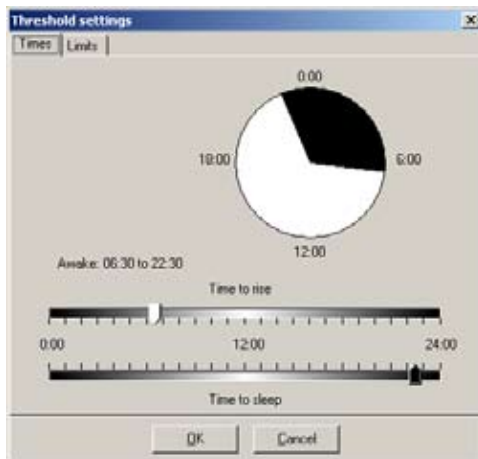


Fig. 10: Awake and Asleep Times

Setting BP Threshold Limits

Users may select the systolic and diastolic thresholds used to calculate the blood pressure load. There are two settings for thresholds: global and custom.

Global thresholds apply to all patient files; custom settings apply to all studies for a particular patient. The default setting for *AccuWin Pro v3* is a global threshold set to 140/90 for awake periods and 120/80 for asleep periods. If desired, you can customize the global thresholds or choose from a published standard: JNC7/AHA or ESH. JNC7¹, or the Seventh Report of Joint National Committee on the Prevention, Detection, Evaluation, and Treatment of High Blood Pressure, published standards are 135/85 for awake periods and 120/75 for asleep periods. The American Heart Association (AHA)² published a normal 24-hour BP standard at 130/80. The ESH³, or European Society of Hypertension, guidelines are 135/85 for awake periods and 120/70 for asleep periods.

1 JNC 7: National High Blood Pressure Education Program. "The Seventh Report of the Joint National Committee on Prevention, Detection, Evaluation, and Treatment of High Blood Pressure". *NIH Publication No. 03-5233* May 2003.

2 AHA: Pickering, TG, et al. (2005). "AHA Scientific Statement: Recommendations for blood pressure measurement in humans and experimental animals". *Hypertension*, 45, 142–161.

3 ESH: O'Brien, et al. (2003). "European Society of Hypertension recommendations for conventional, ambulatory, and home blood pressure measurement". *Journal of Hypertension*, 21, 821–848.

Alternatively, you can define **custom thresholds**, which apply to a particular patient's file. These thresholds will apply to all studies for this patient only.

1. Display an ABP study (Fig. 8). Choose the **Threshold** tab.
2. Click **Edit Thresholds** to change the threshold settings.
3. Click Use custom or Use global thresholds and move the sliders to change the ranges according to your preferences.
4. Click **OK** to apply the changes.

Setting Pediatric Threshold Limits

The American Heart Association recently published a statement supporting ABP monitoring in children, as hypertension is being diagnosed with increasing frequency in pediatric patients.¹ The guidelines for hypertension differ from those for adults; hypertension is generally defined as a blood pressure measurement above the 95th percentile for that gender and height. Specific ABP thresholds are programmed into AccuWin Pro v3.²

To apply a pediatric threshold:

1. Display an ABP study (Fig. 8). Choose the Threshold tab.
2. Click **Pediatric Thresholds** to change the threshold settings.
3. Enter the patient's sex and height. The height may be entered in centimeters (cm), inches, or feet/inches. The height may be entered by typing it directly into the text box or by using the up/down arrows to the right of the text box (Fig. 12).
4. Click **Apply**. Pediatric thresholds are applied as a **custom threshold**, which apply to a particular patient's file.

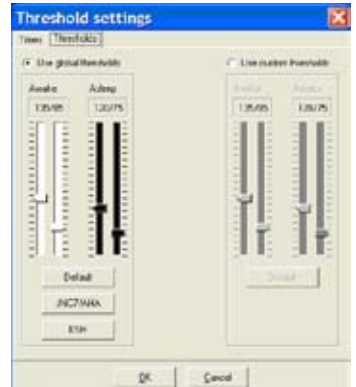


Fig. 11: Setting BP Thresholds

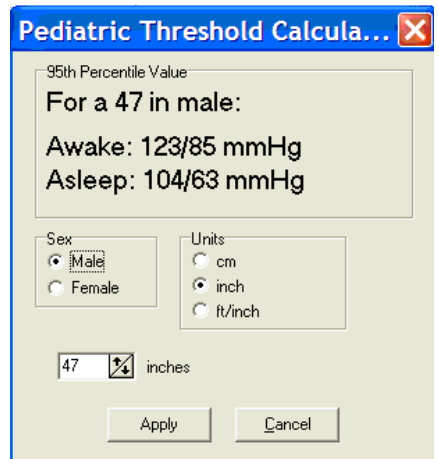


Fig. 12: Pediatric Threshold Calculator

1 Pediatric Threshold Calculator: Urbina, E. et al. (2008). "Ambulatory Blood Pressure Monitoring in Children and Adolescents: Recommendations for Standard Assessment, A Scientific Statement from the American Heart Association, Atherosclerosis, Hypertension, and Obesity in Youth Committee of the Council on Cardiovascular Disease in the Young and the Council for High Blood Pressure Research". *Hypertension*, 52, 433-451.

2 Pediatric Threshold Calculator: Soergel, M. et al. (1997). "Oscillometric twenty-four-hour ambulatory blood pressure values in health children and adolescents: A multicenter trial including 1141 subjects". *Journal of Pediatrics*, 130(2), 178-184.

Defining Time-Slice Periods

Users may define up to twelve specific time windows for which to perform a statistical analysis. The statistical analysis corresponding to time-slice periods will appear in the **Time-Slice Statistics** tab to the right of the **Statistics** tab.

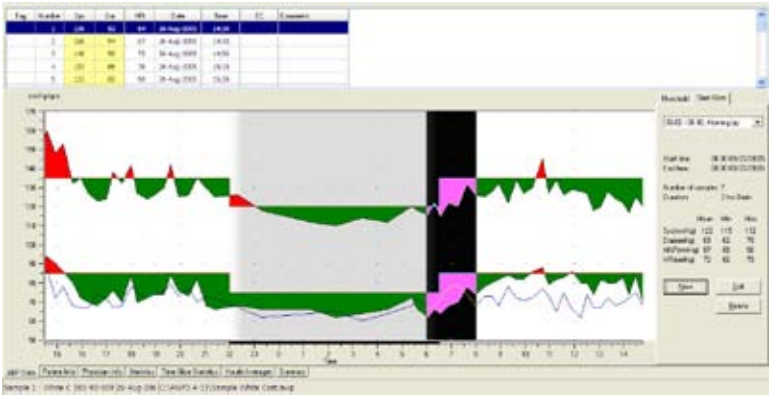


Fig. 13: ABP Data Graph, Time-Slice tab

To create a time-slice:

1. Display the ABP study. Choose the Time-Slice tab (Fig 13).
2. Choose the **New** button. Click on the graph to determine the starting time of the time-slice period, and drag the mouse to the end time. Releasing the mouse button determines the ending time of the time-slice. The time-slice will be highlighted on the ABP data graph.
3. Once the time-slice is created, a **Time-Slice Settings** box (Fig 14) will appear and allow the user to name the time-slice and edit the start and end times. Choose **Save** to create the new time-slice.

To edit a time-slice:

1. Display the ABP study. Choose the Time-Slice tab (Fig. 13).
2. Choose the time-slice from the drop down menu. The time-slice period will be highlighted on the ABP data graph.
3. Choose **Edit**. The **Time-Slice Settings** box (Fig. 14) will appear.
4. Edit the name or start and end times. Choose **Save** to edit the time-slice.



Fig. 14: Time-Slice Settings box

To delete a time-slice:

1. Display the ABP study. Choose the Time-Slice tab (Fig. 13).
2. Choose the time-slice from the drop down menu. The time-slice period will be highlighted on the ABP data graph.
3. Choose **Delete**. Click **Yes** to delete the time-slice.

Editing Patient Info and Physician Info

Patient information and physician information may be entered while data is retrieved from the monitor (see Retrieving Data from the ABP Monitor). Users may edit or add this data once the patient file is saved as well.

To edit **Patient Info**:

1. Choose the Patient Info tab at the bottom of the screen.
2. Click on the **Edit** button.
3. Update any patient information and choose **OK**.

To edit **Physician Info**:

1. Choose the Physician Info tab at the bottom of the screen.
2. Click on the **Edit** button.
3. Update any information and choose **OK**.

Users are able to enter and edit Patient History, Reason for Test, Current Medications, and Physician Interpretation. As an alternative to typing in the free text field for **Current Medications**, users may add **Medication**, **Dosage**, and **Frequency** to dictionaries that are stored within *AccuWin Pro v3*.

To use the **Add Medications** function (Fig. 15):

1. Type a medication into the drop-down titled **Medication**.
2. Once the medication is typed, press "Enter" or "Tab" on the keyboard or use your mouse to move the cursor to the Dosage column.
3. Enter **Dosage** and **Frequency** in the same manner. Pre-populated entries in the **Frequency** column are listed below.

Abbreviation	Meaning
ac	Before meals
bid	Twice daily
pc	After meals
prn	As needed
qam	Every day before noon
qd	Every day
qh	Every hour
qhs	Every night at bedtime
qid	Four times a day
qod	Every other day
qpm	Every day after noon
qwk	Every week
tid	Three times a day

- Once the Medication, Dosage, and Frequency are entered into their respective columns, choose the “Add” button to add the information to the **Current Medications** field above and to the dictionary.

The screenshot shows a form titled "Add Medications". It contains three dropdown menus: "Medication", "Dosage", and "Frequency". To the right of these fields is a rectangular button labeled "Add".

Fig. 15: Physician Information, Add Medications

To locate and choose a previously stored entry:

- Either start to type an entry or choose the down-arrow with your mouse. As more characters of the desired entry are typed, a list of matching entries will be displayed.
- To choose the stored entry, use the down arrow on the keyboard, highlight the desired entry, and press “Enter” on the keyboard. Alternatively, use the mouse to highlight the entry, and then double-click the entry.
- When each of the columns are populated with the correct information, press “Add.”

The screenshot shows a "Patient Information" form with the following fields and values:

Patient name	Sample 1 - White Coat	Birthday	17 Mar 1960
Patient ID	001 00-0000	Sex	Male
Address	12 Parkway Street Chapel Hill, NC 27701 USA	Race	African American
Phone	(919) 444-4543	Height	5' 10"
		Weight	170

At the bottom of the form, there is a tabbed interface with the following tabs: AEP Data, Patient Info, Physician Info, Statistics, Time-Slice Statistics, Hourly Averages, and Summary. The "Physician Info" tab is currently selected.

Fig. 16: Patient Information

AccuWin Pro v3 displays an ABP study in a number of tabular, statistical, and graphical forms that simplify the tracking of a patient's progress and the analysis needed to formulate a diagnosis. Data analysis of the edited ABP study is shown on the following tabs:

- Statistics
- Time-Slice Statistics
- Hourly Averages
- Summary

Viewing Statistics

To view the statistical analysis for the displayed ABP study, click on the Statistics tab on the bottom of the screen. This window shows the following tables:

- **Overall:** Data for all periods combined
- **Awake Period:** Data collected while the patient is awake
- **Asleep Period:** Data collected while the patient is asleep
- **White Coat Period:** If White Coat analysis is enabled on the Summary tab, this table appears and shows data collected during the first hour.

NOTE: All data in these tables is included in exports, except white coat analysis values. White coat analysis values are not included in data exports.

- **BP Load:** Blood pressure load characterizes the distribution of measurements in an ABP study relative to thresholds indicative of hypertension. The percentage of readings above this threshold is the BP load. Formulated to assess the effects of anti-hypertensive therapy in an ABP study, it has been shown to correlate with several indices of hypertensive cardiac target organ involvement.
- **Dip:** The nocturnal fall of blood pressure expressed as a percentage of the mean awake values.



Fig. 17: Figure Statistics

Viewing Time-Slice Statistics

To view the statistical analysis for user defined time-slice periods, click on the **Time-Slice Statistics** tab (Fig. 18).

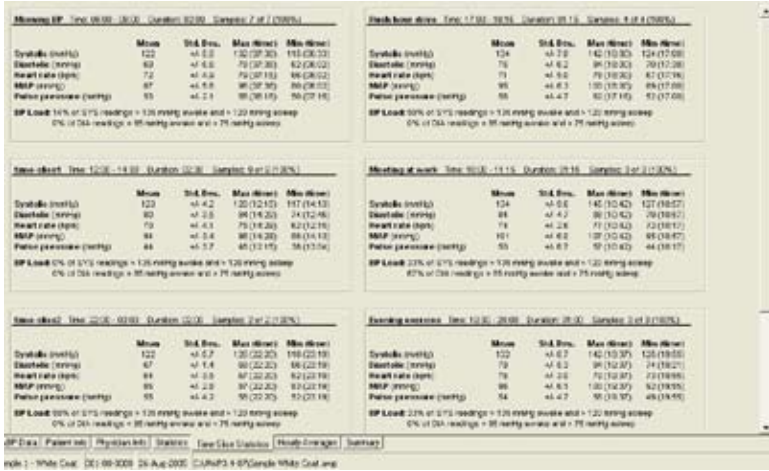


Fig. 18: Time-Slice Statistics

Viewing Hourly Averages

To display trends in a patient's ABP study, click the Hourly Averages tab at the bottom of the window (Fig. 19). The table shows the average BP readings for each hour the patient was tested. In addition to the systolic (Sys), diastolic (Dia) and the heart rate (HR), the following statistical averages are provided:

- **Mean Arterial Pressure (MAP)** - This is the average pressure in an artery over the period of one heartbeat. In the brachial artery, it is calculated by adding the diastolic to one-third of the difference between the systolic and diastolic readings.

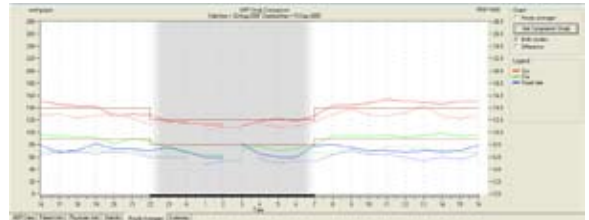


Fig. 19: Hourly Averages

- **Pulse Pressure (PP)** - This is calculated by subtracting the diastolic from the systolic reading. It is another hemodynamic parameter that may serve as an indicator for cardiovascular risk.
- **Pressure Rate Product (PRP)** - This is the product of the average systolic reading multiplied by the average heart rate. PRP strongly correlates to a patient's activity level and may be a key indicator of cardiovascular risk. As typically reported in research and clinical applications, PRP data is divided by 1000.

Use the check boxes at the top of each column to select the statistics you would like to see displayed in the graph at the bottom of the page. You can resize the table and zoom in on the graph with your mouse pointer. (For instructions on resizing and zooming, see Viewing an Ambulatory Blood Pressure Study.)

Viewing the Interpretive Report Summary

To view an interpretive summary for the displayed ABP study, click the Summary tab on the bottom of the screen. This window shows a diagnostic interpretation of the collected ABPM data over a 24-hour period. The report shows awake and asleep systolic/diastolic pressure readings, asleep “dipper” status, and an optional white coat analysis.

Use the Interpretive Summary to evaluate ABPM values across a 24-hour period. Because traditional BP guidelines might not apply to the 24-hour ABPM readings, *AccuWin Pro v3* includes the recommendations for evaluation of ABPM levels according to either:

- Joint National Committee Seventh Report (JNC 7)¹ and American Heart Association (AHA)² guidelines.
- European Society for Hypertension (ESH)³.



Fig. 20: Interpretive Report Summary

AccuWin Pro v3 calculates the asleep dip percentage⁴. Based on the dip percentage, *AccuWin Pro v3* then determines whether the patient is one of the following:

- **Dipper (normal):** A patient shows a decrease of 10% or more systolic or diastolic blood pressure readings taken during sleep compared to readings taken when awake.
- **Non-dipper (abnormal):** A patient shows no decrease or less than a 10% decrease, meaning his or her blood pressure does not fall during sleep compared to the readings while awake.

You can choose to enable or disable the White Coat analysis⁵. The white coat syndrome affects patients whose blood pressure is abnormally elevated when they are in a medical setting.

- 1 JNC 7: National High Blood Pressure Education Program. “The Seventh Report of the Joint National Committee on Prevention, Detection, Evaluation, and Treatment of High Blood Pressure”. *NIH Publication No. 03-5233* May 2003.
- 2 AHA: Pickering, TG, et al. (2005). “AHA Scientific Statement: Recommendations for blood pressure measurement in humans and experimental animals”. *Hypertension*, 45, 142–161.
- 3 ESH: O’Brien, et al. (2003). “European Society of Hypertension recommendations for conventional, ambulatory, and home blood pressure measurement”. *Journal of Hypertension*, 21, 821–848.
- 4 Night-time dip analysis: Claves, H., et al. (2005). “The reproducibility of dipping status: beyond the cutoff points”. *Blood Pressure Monitoring*, 10, 201–205.
- 5 White Coat syndrome analysis: Owens, P, et al. (1999). “Diagnosis of White Coat Hypertension by Ambulatory Blood Pressure Monitoring”. *Hypertension*, 34, 267–272.

Enabling the White Coat analysis option does the following:

- Adds a white coat period to the statistical analysis.
- Displays the first hour of the study, or white coat period, with a light blue background in the graphs.
- Evaluates for and reports on the presence of White Coat syndrome on the interpretive summary, when present.

Comparing two studies

To further enhance the utility of ABP, you may want to measure the progress of a patient's blood pressure by collecting multiple ABP studies. *AccuWin Pro v3* allows comparison of the hourly averages and statistics for two studies.

To Compare Two Studies on the Same Patient:

1. Click on Hourly Averages tab at the bottom of the window.
2. Click on the Get Comparison Study button at the right of the window. (If the button is grayed out, the patient file you are viewing contains only one study.)
3. The Choose Date will appear. Select the study you want to use as a reference for comparison by highlighting it and clicking OK.

The statistics generated during the previous study are now incorporated into the Hourly Averages window. The table expands to include two new lines for each hourly average. One shows the hourly average for the corresponding hour for the reference study, and the other line shows the difference between the displayed and the reference studies.

On the graph at the bottom of the window, both studies' data are displayed (Fig. 21). The dotted lines represent data from the reference study, and solid lines represent data from the displayed study. Select Difference, and the graph will display one line representing the difference between the two studies.

The same navigation tools for the Hourly Average and ABP Data windows also apply to the comparison window. You can move the table to a particular hour by clicking above that hour on the graph. Additionally, you can use your mouse to resize the table and zoom the graph in and out. (For instructions on resizing and zooming, see Viewing an Ambulatory Blood Pressure Study.)

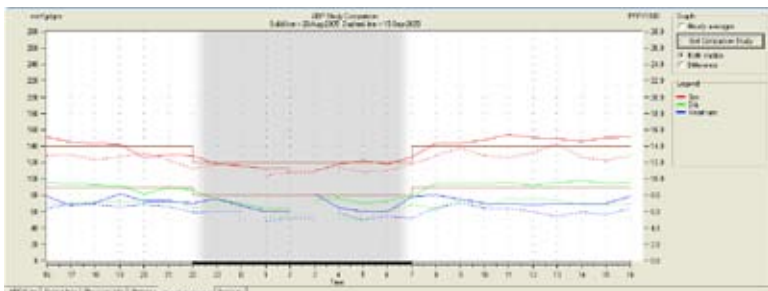


Fig. 21: Hourly Averages, Study Comparison

Creating Reports

Configuring and customizing the report

To document a study and its findings, you can create a customized report. It can contain everything that you have edited and reviewed in the study, or it can comprise a simple, single-page summary report.

AccuWin Pro v3 provides the following preconfigured report formats:

- **One page report:** Includes only the Summary page.
- **Standard report:** Includes the Summary, Patient Information, Statistics, Edited Data, Omitted Data, and Monitor Configuration pages.
- **Full report:** Contains all report pages.

NOTE: All report pages include the *SunTech* logo, patient name, patient ID, and test date.

To configure or customize your report:

1. From the Report menu, select Configure. This displays the Report tab in the Configuration window.
2. On the Report Formats tab, do one of the following:
 - Choose one of the preconfigured formats listed in the left pane. When you select a format, the pages included in that format are selected in the Report Pages box.
 - Choose a custom report format that you created from the left pane.
 - Create a new custom report format by clicking New Report. Enter the report format name. Select the pages from the Report Pages list on the right.
3. To include an additional page in your report, click the page in the Report Pages box. For details, see “Choosing report pages.”
4. Under Report Title, enter text that will print as the header at the top of every report page. For example, enter the physician name and contact information.
5. You may also add a customized logo to appear right-justified in the report header. The logo must be named `report_logo.bmp` and located in the same Windows folder as the *AccuWin Pro v3* executable file. Larger images may produce undesirable effects.

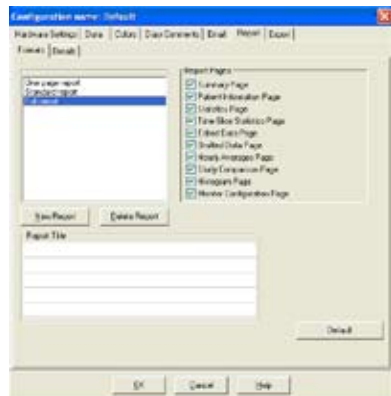


Fig. 22: Configuration, Report Formats



Fig. 23: Configuration, Report Details

To configure report details:

1. Click the Details tab to choose
 - **Graph Type:** Select (for printed report) line, bar, or the type of graph currently displayed.
 - **Graph Scale:** Select 12, 24, or 48-hour time scale for the Summary Report ABP graphs.
 - **Summary Page User Fields:** The page includes two additional sections of your choice. Select one of the following options that includes the two sections:
 - **Medications/Physician Interpretation:** The summary page will include sections for the patient's medications and the physician's interpretation.
 - **Patient History/Physician Interpretation:** The summary page will include sections for the patient history and the physician's interpretation.
 - **Test Reason/Physician Interpretation:** The summary page will include sections for the reason for the test and the physician's interpretation.
 - **Medications/Patient History:** The summary page will include the patient's medications and the patient's history.
 - **Medications/Test Reason:** The summary page will include the patient's medications and the reason for the test.
 - **Patient History/Test Reason:** The summary page will include the patient's history and the reason for the test.
 - **Horizontal Sight Lines:** Select frequency of sight lines on Edited Data or Omitted Data pages.
2. Click OK to apply the changes.

Report - Formats tab

You can choose from the following pages to customize your report. All pages include a header identifying the patient name, ID, and test date. Customize the header with the address of your clinic or any information that you would like to display; this shows as your report title. The footer of each page identifies the type of monitor used and the page number.

The report prints out in the following order:

- **Summary:** The Interpretive Summary appears above the data graph and summary statistics. Summary statistics for the overall, awake, asleep, and white coat (if chosen) periods include: mean systolic and diastolic measurements, mean heart rate, bp load, and asleep dip. Two sections, such as Patient History and Physician Interpretation, selected from options on the Report Configuration Details tab are also included. The referring physician and interpreting physician names, signature lines, and date lines appear at the bottom of the Summary page.
- **Patient Information:** Includes patient name, ID, address, and physical

description. Also displays the two user-selectable fields that were not chosen on the Summary page. These could include patient history, reason for test, current medications, and physician interpretation.

- **Statistics:** Displays detailed statistics, including BP loads, averages, standard deviations, minimums (with the time of occurrence in parentheses), maximums (with the time of occurrence in parentheses), and pulse pressure for each period and for the entire study. The table displays data by different time periods: overall, awake, and asleep. The Overall section includes the dip percentage. Data collected during the White Coat period (the first hour of collection) appears if the White Coat Analysis option was enabled on the Summary tab.
- **Time-Slice Statistics:** Displays detailed statistics, including BP loads, averages, standard deviations, minimums (with the time of occurrence in parentheses), maximums (with the time of occurrence in parentheses), and pulse pressure for each user defined time-slice period.
- **Edited Data:** Shows all the BP readings used in calculating statistics in graph and table format. Displays the time of the reading, the BP measurement, heart rate, mean arterial pressure (MAP), codes, and comments. The vertical black line between # and Time columns indicate the sleep period. A ">" symbol adjacent to the BP reading indicates that it has exceeded the threshold.
- **Omitted Data:** Shows the BP readings that were not used to calculate statistics. This page appends to the Edited Data page if printed.
- **Hourly Averages:** Displays hourly average data for the study in graph and table formats. In addition to the BP reading, heart rate, and MAP, Pulse Pressure (PP) and Pressure Rate Product (PRP/1000) are included. Prints as displayed on the Hourly Averages tab.
- **Study Comparison:** Displays study comparison data in graph and table formats. Prints as displayed on the Hourly Averages tab.
- **Histogram:** Contains frequency histograms for all variables.
- **Monitor Configuration:** Shows the study parameters used by the ABP monitor for the study. The page shows information about sleep times and definitions for applicable codes. This page also displays the selected maximum pressure value.

Report - Details tab

You can control the following parameters for printing reports:

- **Graph Type:** Controls if the type of graph printed will be line, bar, or the graph type currently displayed.
- **Graph Scale:** Controls the time scale for the graph of ABP data when printed. This graph will be printed when the ABP summary page, Edited Data Page, or Omitted Data Page are printed.
- **Horizontal Site Lines:** Sets how often sight lines are used when printing the Edited Data Page, or Omitted Data Page in the report.

Previewing the report

1. To view individual pages after the report has been configured, from the Report menu, select Preview. Alternatively, click the Print Preview toolbar button. The Print Preview window is displayed (Fig. 24).
2. To view individual pages, click on the numbered icons on the left side of the window.
3. To change the format and details, click Configure.
4. To print the report, click Print Report. To print the currently viewed page, click Print Page.

Printing the report

AccuWin Pro v3 provides you with flexibility in printing the report. You can either:

- Print the configured report by clicking the Print button on the toolbar or by going to the Report menu and selecting Print. The Print dialog box is displayed, where you can choose the printer under Name.
- Print the report page corresponding to the current tab that is being displayed by clicking Print Page. The table below shows the report page(s) generated by the display of each tab.

Displayed Tab	Report Page(s) Printed
ABP Data	Edited and Omitted Data with Graph
Patient Info	Patient Information
Physician Info	Summary Page
Statistics	Statistical Tables
Time-Slice Statistics	Statistical Tables
Hourly Averages	Hourly Averages and Study Comparison
Summary	Summary Page

Saving the report as a PDF

To save the report as an Acrobat Portable Document Format (PDF), from the Report menu, select Create PDF. Alternatively, click the PDF toolbar button. The system creates the PDF to the location you specify.

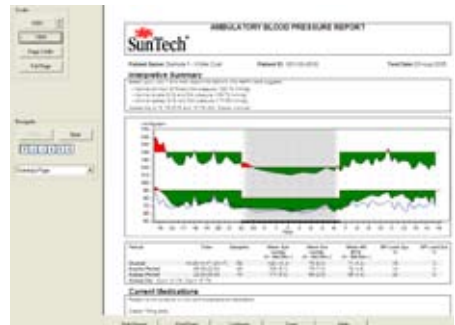


Fig. 24: Print Preview

NOTE: You must have a printer installed. If you do not have a printer installed, see the online Help to install a driver to print to a file.

Managing Patient Files

Manage patient files using the following functions.

Open

This feature allows you to retrieve a patient's data file.

1. From the toolbar, select the BP Data button or from the File menu, select Open. The Open APB Study – Choose Patient dialog box appears, allowing you to select the patient file to open.
2. Select the patient whose BP data you wish to retrieve by clicking on the patient's name. Select OK.
3. If the patient has more than one ABP study, a second dialog box, Choose Date, will be displayed. Select a date and press OK.
4. The patient's data appears in the Display Area. You can now edit the BP data or print a report.

NOTE: If patient files are saved to a shared folder/network, the patient file may be accessed by only one user at a time. If an additional user attempts to open a patient file in use, a message stating that the file is currently in use will be displayed.

Copy

To copy BP data files to other disks or directories:

1. From the File menu, select Copy.
2. The Copy – Choose File(s) dialog box appears. Select the files you want to copy. Click OK.
3. A Copy Files to: dialog box appears. Choose the destination drive and directory. Click OK to copy the files.

Export

To export BP data files to either ASCII or GDT format:

1. From the File menu, select Export.
2. Choose either to export the Displayed File or select files From Database.
3. If selecting files from the database, check the files you want to export and click OK.
4. The Export dialog box appears. Click Export to export the file(s). Select Edit to change the files you want to export. Choose Configure to configure the details of the exported file(s). (For instructions on configuring export, see Export in Customizing and Configuring *AccuWin Pro v3*)

The exported files will be named after the *AccuWin Pro v3* filename as either a *.ASC or *.GDT file depending on the format chosen.

Email

To email *AccuWin Pro v3* data files as an attachment:

1. From the File menu, select Email.
2. Choose either to email the Displayed File or select files From Database. If Displayed File is chosen, the data currently displayed on the screen is emailed.

If From Database is chosen, the patient database window opens and prompts you to Select file(s) to email. All selected tests will be emailed.
3. If selecting files from the database, check the files you want to email and click OK.
4. The Email dialog box will appear. Type in the recipient email address(es), a subject for your message, and a message. The *AccuWin Pro v3* file must have a ".awp.enc" file appendage for the recipient to view the *AccuWin Pro v3* file. The ".enc" denotes that the file is encrypted.
5. Click Send to email the file(s). Select Attach to change the files you want to email. Choose Configure to configure the details of emailing the file(s). (For instructions on configuring email, see Email in Customizing and Configuring *AccuWin Pro v3*).

Delete

This feature allows you to delete one or more patient data files. Each patient file contains ALL the BP data for that patient; therefore ALL BP data for that patient will be deleted.

1. From the File menu, select Delete.
2. A Delete – Choose File(s) dialog box appears. Select which BP data files to delete. Select or de-select items by clicking in the box beside the patient name. After you select all the files to be deleted, click OK to continue.
3. A dialog box to confirm the deletion of these files is displayed. Select Yes to continue.

Converting patient data files

If you have patient files from previous versions of *AccuWin Pro v3* and you try to open or copy them, a message appears requesting that you convert them to new encrypted files to enhance security of patient information.

To convert patient data files:

1. Begin to open or copy a patient file saved in a previous version of *AccuWin Pro v3*. A message appears indicating that one or more patient files is not in a HIPPA compliant format and requests that you convert them now.
2. Do one of the following:
 - To convert all patient files at this point, we recommend that you first back up these files. Then, click Yes.
 - To not convert all patient files now, click No. If you choose not to convert them, an asterisk (*) appears next to the patient file name in the list of files to indicate files that have not been converted as encrypted files. Each time you launch *AccuWin Pro v3* and open patient files saved in the previous format, the message to convert the files appears.

Customizing and Configuring *AccuWin Pro v3*

In one facility it might be helpful to have multiple configuration profiles for your users, for example, one for each physician in the practice or one for each service bureau client.

You can create or edit configuration profiles in *AccuWin Pro v3*:

1. From the Configure menu, select Preferences. Create a new profile or edit an existing profile.
2. Click on the name of the profile you want to edit and click Edit.
3. To create additional profiles, click New, type in a name to identify the profile, click Enter, and then select Edit.
4. Click OK to choose a configuration

If you are editing a configuration, a window opens that displays a number of functions shown as tabs across the top of the window. These are described in detail on the following pages.

NOTE: On each tab, clicking a Default button returns the settings of that tab to the default factory settings.

Hardware Settings

This window (Fig. 26) allows you to configure your computer to communicate with your ABP monitor. You may also test the settings that you have selected to make sure it is operational. For complete instructions, see *Configuring your computer for Communication*.

1. Under ABP device, select the unit.
2. Under Serial port, select the COM port into which you have plugged the unit.
NOTE: This is still required if using the USB adapter.
3. Click Test port. The unit will beep once, and you will receive a verification message that the port is properly configured.
4. If you receive an error message, "Cannot communicate with ABP device," one of the following applies:
 - You have not selected the appropriate port or device.
 - The device is not connected properly.
 - There is no battery power in the monitor.
 - The USB adapter has not been installed correctly.

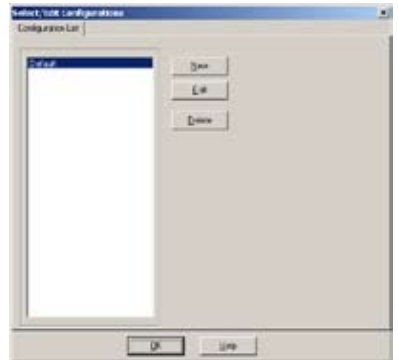


Fig. 25: Configuration list



Fig. 26: Configuration, Hardware

Data Settings

This window (Fig. 27) includes the following features.

- **Default path:** This is the folder created on your hard drive when you installed *AccuWin Pro v3*. If you do not want to save reports in this folder, select another in the tree (picture below) and it will automatically become the default. If you are creating multiple configurations, you may want to set up a separate folder for each profile. You can choose to have *AccuWin Pro v3* always access this folder, access it on the first search only, or access the last folder used. Change this by selecting one of the three choices shown.
- **Always prompt for study start time/date:** Check this to have the software prompt you to enter the test date when you retrieve data collected by an ABP monitoring device.
- **Show hints:** Check this to view information about icons and graph features as you work.
- **Auto-detect serial ports:** Check this to have the system check for serial port connections.
- **Enable unique file identifier:** Check this to add a unique file prefix to the ABP-assigned file name. This option adds a user-defined four-character alphanumeric prefix to the file name (for example, xxxx-yyyy-mm-dd-##.awp where xxxx is the user defined prefix.)
- **Enable prompt for file name:** Check this to have the system prompt you for a file name when you save the ABP file.

NOTE: If both options are selected, the prefix is added to the user-entered file name (for example, xxxx-<any user-defined prompt>).

To return all settings in this window to the factory defaults, click Default.

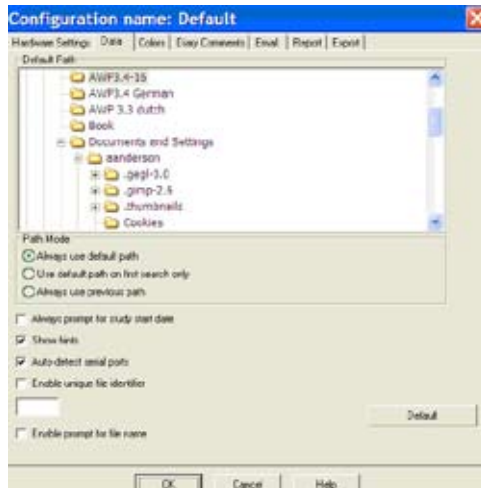


Fig. 27: Configuration, Data

Color Settings

This feature allows you to customize the color scheme of the graphs and charts generated through each configuration profile you have created (Fig. 28). In large organizations, varying the color palette for each profile allows printed reports and onscreen data to be easily matched to the correct physician or client. Customizable palettes also allow you to fine-tune the color scheme for optimum reproduction by your office printers. Shading on the ABP graph indicates the period of the study. The default colors (which can be changed using the Configuration > Colors tab) are the following:

- **White coat period:** Shown in light blue shading and indicates the white coat period, which is the first hour of the study.
- **Sleep period:** Shown in light gray shading and indicates the asleep period of the study.

Here are the steps for customizing the color palette for a configuration profile:

1. Click on the tab for the attribute you want to change. In the example shown here, the ABP graph tab is selected, but this function works identically for each tab.
2. To change a color, click on it. A color palette window will appear, showing 48 basic colors. You may select one of these, or you may create your own custom colors by clicking on the Define Custom Colors button.
3. Repeat this process for each color you want to change.

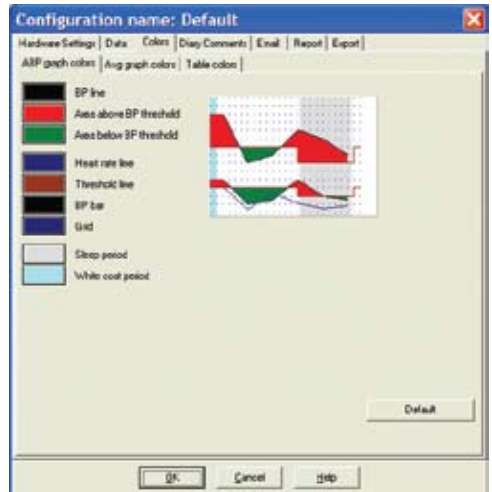


Fig. 28: Configuration, Colors

Diary Comment Settings

You can use this feature to keep track of entries from patient diaries describing the patient's activities during BP measurements. The entries on the default list in this window (Fig. 29) correspond to the lists contained in patient diaries. You can customize this list by adding new entries or deleting some of the default entries.

The entries on your final list are incorporated into a list that appears in the Comments column on the table located on the ABP Data tab. When you are processing a file, you can select one of the entries from this menu with your mouse instead of typing the text. (See Patient Diary Comments.)

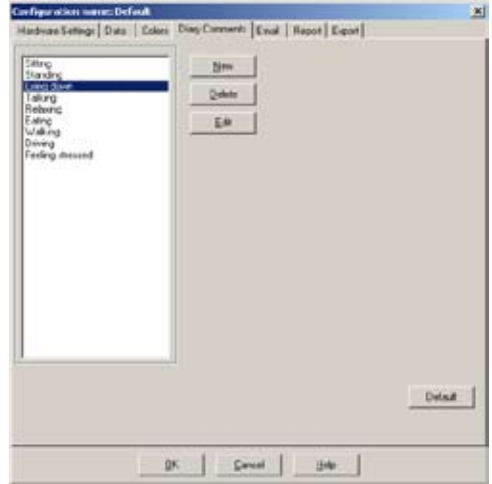


Fig. 29: Configuration, Diary Comments

Email Settings

The Email feature creates a new email message, attaches selected patient data files, and sends the email to its destination.

To set up the parameters for email (Fig. 30):

1. Choose the type of connection you are using by selecting LAN or Dial up.
2. Fill in the Email server settings. This information is available from your Internet service provider (ISP) or office network administrator. Your ISP may require you to check the Authorize box in order for you to send email.
3. The Email settings you fill in will appear in the header of your email: the name of the recipient, your return email address, the subject of the email and a message of your choosing.

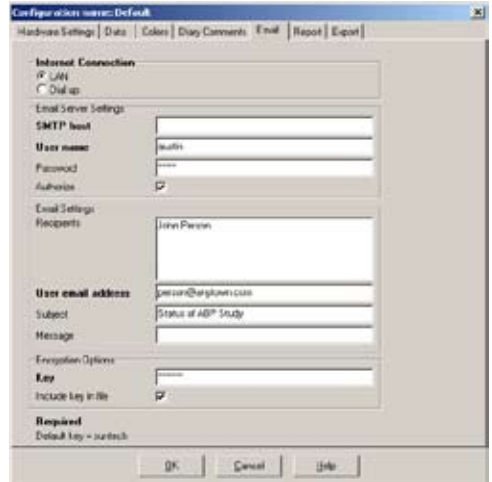


Fig. 30: Configuration, Email

Encryption Options allow you to encrypt your emails with a key of your choosing, which you will need to make available to the recipient of the email. Patient files are always encrypted. If the key is included in the file, any copy of AccuWin Pro v3 can open the file. To increase security, uncheck the box for Include key in file and type in an encryption key. The recipient of the emailed file must have the key to open the attached file.

All bolded categories on the tab are required.

Report Settings

This feature allows you to format your reports and to choose what information is included in them.

To format a report:

1. From the Configure menu, select Preferences. Select the profile you want to edit and click Edit.
2. Click the Reports tab and then the Formats tab.
3. Select which items will be included in your One Page report, Full Report, and Standard Report. Select or de-select the check-boxes on the right to make these changes.
4. To create a new report, click New Report, enter a name of your choosing, and select the items you want to appear in that report.
5. To delete any of these reports, select the report and click Delete Report.
6. Enter a Report Title by typing the title of your choice in the box provided. The title appears as a header on each page of the report. Each report format allows for a different report title.
7. As with the other features, return to factory settings by clicking Default.
8. To change the graph type, graph scale, and report summary for your report, click Details. You can also change the number of horizontal sight lines shown in your report to make it easier to read, or you can choose not to use them.

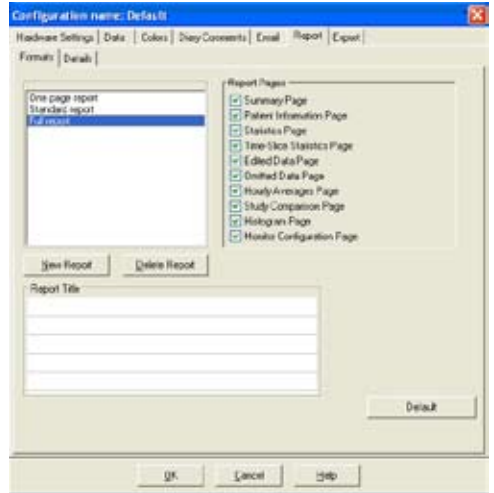


Fig. 31: Configuration, Report Formats

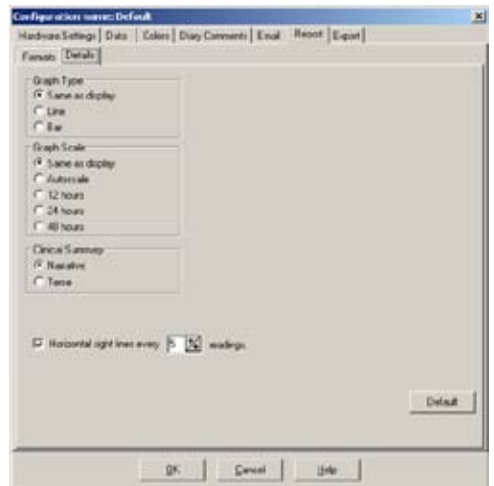


Fig. 32: Configuration, Report Details

Export Settings

This feature allows you to export reports in an alternate format (Fig. 33), which is useful when the data will be merged with a database. You can export reports in either ASCII or GDT output. GDT is a format used to share computer medical records.

1. To configure exports to ASCII output, select options in these areas:
 - **Delimiters:** Choose which character will separate field items in the database.
 - **Export:** Choose exactly what information will be exported.
2. To configure this function for GDT output, do the following:
 - Select the Export GDT tab.
 - Click the GDT Summary checkbox.

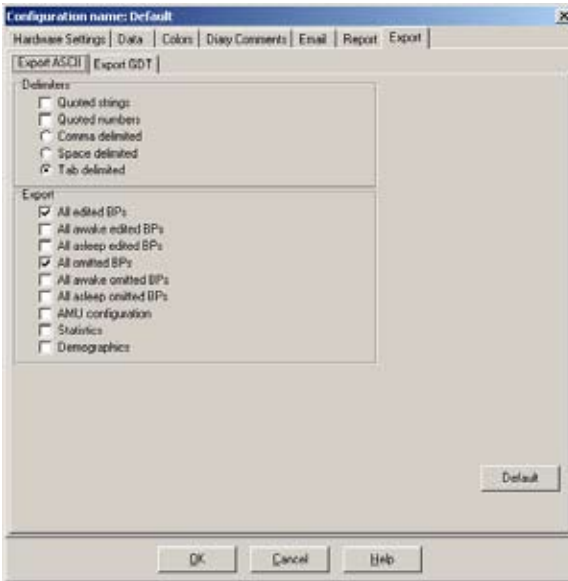


Fig. 33: Configuration, Export

Overview of AccuWin Pro v3 Report Symbols

Report Symbols

The following symbols are used in the Edited and Omitted BP data report pages.

Symbol	Explanation
"†"	Indicates a reading was taken in the Office Test Mode.
"+"	Indicates that the reading was initiated by pressing the Start/Stop key on the monitor.
"-"	Indicates that the pressure to which the cuff inflated in the dynamic inflation mode may have been too low to obtain an accurate systolic reading. Readings marked with "-" are automatically omitted from the <i>AccuWin Pro v3</i> report. Compare the "tagged" BP reading marked with a "-" to the BP readings above and below it. If the omitted reading seems to be reasonable, "un-tag" it; if it does not seem reasonable, leave it "tagged".
"r"	Indicates that this reading was a retry reading. Retries are readings that are automatically obtained 4 minutes after a questionable reading is obtained by the monitor.
"!"	The exclamation point is a permanent tag and is used when an accurate reading can not be determined by the ABP monitor.
">"	This symbol will appear to the left of the measurement if the reading was above the set threshold.

Overview of Event Codes

Event (EC) codes are used extensively in your decisions on editing BP data. The codes are displayed on the tables in the ABP Data tab and in the Edited and Omitted BP data report pages in the columns labeled EC.

- Event Codes describe the conditions under which the BP measurement was taken. They are also accompanied by solutions that can be taken to avoid their occurrence in the future.

NOTE: Codes do not automatically indicate that a reading is invalid; they serve only as guides to help you edit data.

Event Code Definitions

Event Code	Description in AccuWin Pro v3	Solution
1	Weak or no oscillometric signal	Check the position of the cuff; tighten the cuff.
2	Artifact/erratic oscillometric signal	Remain still during the BP reading.
3	Exceeded retry count (4 attempts)	Remain still during the BP reading.
4	Exceeded measurement time limit (120 seconds)	Check the air hose connections and make certain the cuff is tight.
85	Reading aborted (blocked valves or pneumatics)	Check the air hose connections and make certain the air tubing is not crimped.
86	Reading aborted (user)	Push the Start/Stop button to restart the reading.
87	Reading aborted (inflate timeout or air leak)	Check the air hose and cuff.
88	Reading aborted (safety timeout)	Retry the reading by pushing the Start/Stop button. If the problem persists, return the unit for service.
89	Reading aborted (cuff over-pressure)	Check for blocks or kinks in the air hose.
90	Service required (power supply out of range or other hardware problem)	Replace the batteries. If the problem persists, return the unit for service.
91	Service required (safety override fitted or autozero out of range)	Retry the reading by pushing the Start/Stop button. If the problem persists, return the unit for service.
97	Service required (transducer out of range)	Return the unit for service.
98	Service required (A/D out of range)	Return the unit for service.
99	Service required (EEPROM calibration data CRC failure)	The unit needs to be recalibrated. Return for service.

Administrative Tools

AccuWin Pro v3 includes administrative tools that implement additional security-related features. With *AccuWin Pro v3* Admin Tools, you can create new user accounts, change user access levels (administrative or user access), reset user passwords, enable the login process, and enable automatic logging off after a specified time elapses without user activity. Access to Admin Tools requires administrative level login.



Fig. 34: Admin Tools window
(admin user name removed)

Maintaining login options

AccuWin Pro v3 includes the flexibility of enabling or disabling password login security. Enabling this option supports the HIPAA access protection requirements. If your facility uses an Electronic Medical Record (EMR) system or another system that already provides login security, you may not want to enable the password control feature in *AccuWin Pro v3*.

When the password control feature is enabled, you can also enable an automatic logoff feature that logs out a user who has not performed any activity in *AccuWin Pro v3* after a specified time has elapsed. You can enable or disable automatic logoff and you can set the time.

NOTE: All Admin Tools options require access from users who have administrative privileges.

To enable login access:

1. From the Configure menu, select Admin Tools. If you have not logged in as an administrator or if the password control feature is not enabled, the Enter Name and Password dialog window appears. Enter your administrative name and password and click OK. The Admin Tools window (Fig. 34) appears.

NOTE: If this is the first time accessing Admin Tools, you must use the default user name and password found on your *AccuWin Pro v3* install CD (Fig. 35).

2. Click the Use Passwords checkbox. A check indicates that the option is enabled.
3. Click Done.

NOTE: You must create a new Administrative user name and password before deleting the default administrator.

To enable automatic logoff:

NOTE: The password control feature must be enabled.

1. From the Admin Tools window, click the Auto Logoff Enabled checkbox. A check indicates that it is enabled.
2. In the Auto Logoff Minutes box that appears after you check the Auto Logoff Enabled checkbox, enter the number of inactive minutes of inactivity that must elapse before the user is automatically logged off. Alternatively, you can click the Up or Down arrows to select a number. You can set the number of minutes between 1 and 240 minutes; the default is 20 minutes.
3. Click Done.

AccuWin Pro v3 includes security features that allow only an administrative user maintain user account information. The administrator can add, modify, and delete user accounts.

There are two types of user accounts: administrator and user. The administrator has access to administrative features such as enabling login access and enabling the automatic logoff feature.

The administrator can also reset user passwords for individual users who have forgotten their passwords.



Fig. 35: AccuWin Pro v3 CD with User Name and Password location. Check your CD label for actual default user name and password.

NOTE: All Admin Tools options require access from users who have administrative privileges.

Adding user accounts

You can add new user accounts quickly and easily. An administrator can create a user at either an administrative or user level. Each new user is assigned a default password of "user". When new users first log in, they must change the default password before they can perform any tasks in *AccuWin Pro v3* (Fig 36). This applies to both administrative and user accounts.

1. From the Configure menu, select Admin Tools.
2. If you have not logged in as an administrator or if the password control feature is not enabled, the Enter Name and Password dialog window appears. Enter your administrative name and password and click OK. The Admin Tools window appears.
3. Click New User.
4. In the User Name field that appears, enter the new user login name. The newly added user automatically has his or her password set to the default "user". When logging in for the first time, the new user will be prompted to change the password.
5. Click OK.
6. To change the level of user, select the user then click Change Level. This option toggles the user level that appears in the Level column between User and Admin.
7. Click Done.



Fig. 36: Admin Tools window, adding user

Changing user level

You can change user account levels to either administrative or user privileges. You can not change user names.

1. From the Configure menu, select Admin Tools.
2. If you have not logged in as an administrator or if the password control feature is not enabled, the Enter Name and Password dialog window appears. Enter your administrative name and password and click OK. The Admin Tools window appears.
3. Click the user account you want to change.
4. To change the level of user, click Change Level. This option toggles the user level that appears in the Level column between User and Admin.
5. Click Done.

Deleting user accounts

You can delete both administrative and user accounts; however, you can not delete an administrator account if it is the only existing one. You also cannot delete yourself. Deleting a user account prevents the user from logging into *AccuWin Pro v3*, if the password control feature is enabled.

1. From the Configure menu, select Admin Tools.
2. If you have not logged in as an administrator or if the password control feature is not enabled, the Enter Name and Password dialog window appears. Enter your administrative name and password and click OK. The Admin Tools window appears.
3. Click the user account you want to delete.
4. Click Delete User.
5. A confirmation message appears. Click Yes to continue deleting the user account.
6. Click Done.

Resetting user passwords

When a user forgets his or her password, that user must contact the system administrator, who will reset the user's password to the default password of "user." When he or she first logs in again, the user must change the default password before being able to perform any tasks in *AccuWin Pro v3*.

1. From the Configure menu, select Admin Tools.
2. If you have not logged in as an administrator or if the password control feature is not enabled, the Enter Name and Password dialog window appears. Enter your administrative name and password and click OK. The Admin Tools window appears.
3. Click the user account you want to change.
4. Click Reset Password.
5. Click Yes to confirm resetting the password.
6. A confirmation message appears, click OK.
7. Click Done.

AccuWin Pro v3 Technical Assistance

Windows 2000 Installation

1. Install AWP using the installer provided on the CD.
2. Open the CD directory and place a copy of the ftd2xx.dll file into the AWPv3 directory (C:\AWPv3.4).
3. Plug in the SunTechMed USB cable to “install” the hardware.
4. The application can now be used.

For any questions, please reference the Help section within *AccuWin Pro*, the product tutorial on the *AccuWin Pro v3* CD, or our website for constant product updates. The FAQ's section on the *SunTech Medical* website can be found here: (http://www.suntechmed.com/bp-devices-and-cuffs/ambulatory-blood-pressure-monitoring/oscar-2#Product_FAQs). If these do not fully address your problem, please contact our service department.

For more detailed assistance, please refer to this *AccuWin Pro v3* User's Guide or contact your nearest *SunTech* service center.

USA, Canadian, and Latin American Headquarters and Service

SunTech Medical, Inc.
507 Airport Blvd, #117
Morrisville, NC 27560-8200

Toll free USA and Canada: 800.421.8626
Phone: 919.654.2300
www.SunTechMed.com

European and Far East Headquarters and Service

SunTech Medical, Ltd.
Oakfield Industrial Estate
Stanton Harcourt Road
Eynsham, Oxfordshire OX29 4TS
+44 (0)1.865.884.234

Asia

SunTech Medical, Ltd.
Level 30,
Bank of China Tower
1, Garden Road Central
Hong Kong
+852.2251.1949

Safety

SunTech Medical™ is responsible for the safety, reliability, and performance of its equipment provided that the following conditions exist:

Assembly operations, extensions, readjustments, modifications and repairs are carried out by *SunTech* authorized personnel. The electrical installation of the relevant room complies with all Local, State, and Federal building and electrical codes. All equipment is used in accordance with instructions for its use.

Contact *SunTech Medical* for information before connecting any devices not mentioned in this *AccuWin Pro v3 User's Guide*.

To ensure patient safety, use only parts and accessories manufactured or recommended by *SunTech Medical*.

A**ABP Study**

Awake/Asleep times	18
Editing	16
Opening patient file	14
Programming	12
Retrieve Data	13

AccuWin Pro

Color Settings	36
Data Settings	35
Diary Comments Settings	37
Email Settings	37
Export Settings	39
Menu Options	10
Configure	11
File	10
Help	11
Monitor	10
Report	11
Report Settings	38
Report Symbols	40
Toolbar	9

Administrative Tools

Adding accounts	45
Default user name	44
Deleting accounts	46
Login options	43
Reset Passwords	46
User level	45

C**Configuration**

Connecting Device	7
-------------------	---

Creating Reports

Configuration	27
Details	29
Formats	28
Previewing	30
Printing	30
Saving in PDF	30

D

Diary Comments	17, 37
----------------	--------

E

Event Codes	42
-------------	----

H

Hardware Requirements	6
-----------------------	---

Settings	34
----------	----

I**Identification**

Login	8
Password	8

Introductions

AccuWin Pro	5
Ambulatory BP Monitoring	4

M

Medication Function	21
---------------------	----

P**Patient Info** 21

Converting data files	33
-----------------------	----

Managing	31
----------	----

Physician Info 21**Q**

Quality Codes	41
---------------	----

R**Report Summary** 25

Comparative Analysis	26
----------------------	----

Symbols	40
---------	----

S**Safety** 48**Service Centers** 47**Software**

Components	6
------------	---

Installation	6
--------------	---

Requirements	6
--------------	---

Statistics

Hourly Averages	24
-----------------	----

Overall	23
---------	----

Time-Slice	24
------------	----

T**Technical Assistance** 47**Test Codes** 41**Threshold Limits**

BP Threshold	18
--------------	----

Pediatric Threshold	19
---------------------	----

Time-Slice Periods 20**U**

USB Connection	7
----------------	---



SunTech Medical, Inc.
507 Airport Boulevard, Suite 117
Morrisville, NC 27560-8200 USA
Tel: + 1.919.654.2300
1.800.421.8626
Fax: + 1.919.654.2301

SunTech Medical, Ltd.
Oakfield Estate
Eynsham, Oxfordshire
England
Tel: +44 (0) 1865.884.234
Fax: +44 (0) 1865.884.235

SunTech Medical, Ltd.
Level 19, Two, International Finance Centre
8, Finance Street, Central
Hong Kong
Tel: +852.2251.1949
Fax: +852.2251.1950

