


Vitals Management Software
for the 247 BP + Vitals

 SunTech® *VitalShare v1*TM
VITALS MANAGEMENT SOFTWARE

Digitize Your Spot Check Vital Signs Workflow

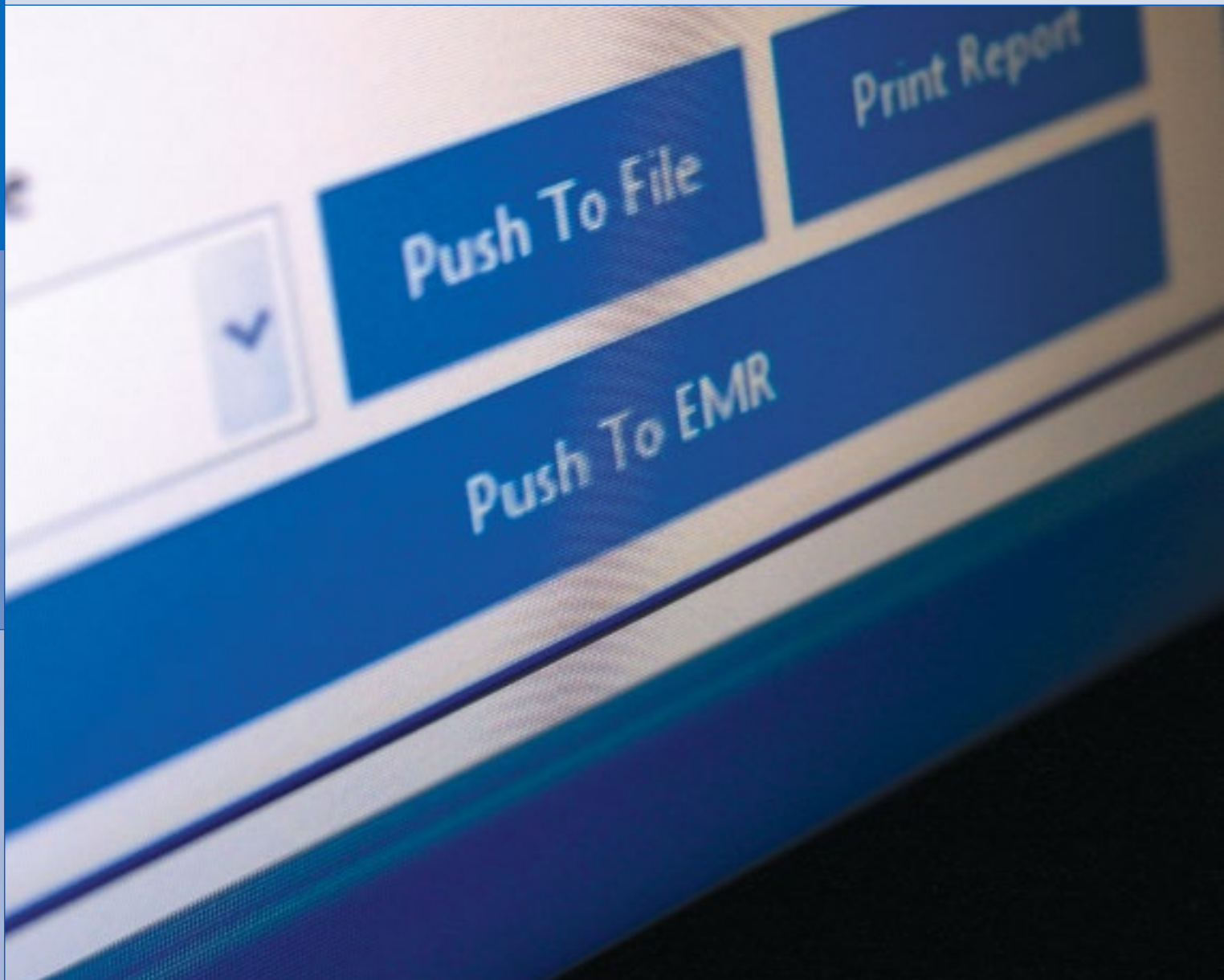


- + Efficient transmission of data with only two clicks
- + Included with all Bluetooth 247 systems
- + Implements IHE Technical Framework profiles using HL7 messaging
- + Supports the validation of patient data

SunTech®'s VitalShare™ software provides a seamless interface between the 247™ BP+Vitals with Bluetooth® and a patient's electronic medical record. *VitalShare* enables wireless communication between the 247 and a networked PC. Clinicians can match vital signs measurements with specific patients and send the data to an EMR without having to write anything down.

+ System Requirements

- PC with *Bluetooth* connectivity
- Windows XP, Windows Vista & Windows 7
- 100 MB free disk space.
- TCP/IP network connection between PC and EMR

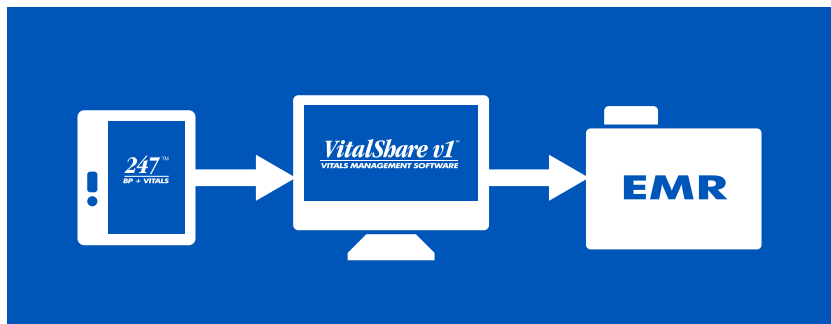

+ IHE International, Inc.

IHE International, Inc. is a multi-stakeholder effort started in 1997 to promote the effective use of electronic health records and health information technology (HIT) by driving the adoption of standards-based solutions for interoperability. IHE's standards-based profiles make sharing information between devices and EMR systems easier.

SunTech's *VitalShare* vitals management software implements all transactions required in the IHE Technical Framework to support the IHE Integration Profiles, Actors and Options listed below.

Integration Profile	Actor	Integration Profile Option
Consistent Time	Time Client	None
Device Enterprise Communication	Device Observation Reporter	None
Patient Demographics Query	Patient Demographics Consumer	None

+ Part Numbers

Description	SpO ₂ Opt.	PN (<i>Blue-tooth</i> ®)
Battery Mobile System, BP, Temperature & SpO ₂	vSat	99-0116-00
	Nellcor®-comp.	99-0089-00
Battery Mobile System, BP		99-0090-00
Battery Mobile System, BP & Temp.		99-0091-00
Battery Mobile System, BP & SpO ₂	vSat	99-0117-00
	Nellcor®-comp.	99-0092-00
Wall System, BP, Temperature & SpO ₂	vSat	99-0118-00
	Nellcor®-comp.	99-0093-00
Wall System, BP		99-0094-00
Wall System, BP & Temperature		99-0095-00
Wall System, BP & SpO ₂	vSat	99-0119-00
	Nellcor®-comp.	99-0096-00
Battery Tabletop System, BP, Temperature & SpO ₂	vSat	99-0120-00
	Nellcor®-comp.	99-0097-00
Battery Tabletop System, BP		99-0098-00
Battery Tabletop System, BP & Temp.		99-0099-00
Battery Tabletop System, BP & SpO ₂	vSat	99-0121-00
	Nellcor®-comp.	99-0100-00



*All Bluetooth equipped systems include a copy of *VitalShare* vitals management software.

Nellcor® is a registered trademark of Nellcor Puritan Bennett LLC

The *Bluetooth* word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by SunTech Medical, Inc. is under license.

SunTech Medical, Inc.
507 Airport Boulevard, Suite 117
Morrisville, NC 27560-8200 USA
Tel: + 1.919.654.2300
1.800.421.8626
Fax: + 1.919.654.2301

www.SunTechMed.com

SunTech Medical, Ltd.
Oakfield Estate
Eynsham, Oxfordshire
England
Tel: + 44 (0) 1865.884.234
Fax: + 44 (0) 1865.884.235

www.SunTechMed.com

SunTech Medical, (Shenzhen) Co., Ltd.
2/F of Building A, Jinxiongda Technology Park
Guanlan, Bao'an District Shenzhen, 518110, PRC
Tel: + 86-755-2958 8810
+ 86-755-29588986 (Sales)
+ 86-755-29588665 (Service)
Fax: + 86-755-2958 8829

www.CN.SunTechMed.com